

VTM-N Collection and Transport Kit

On the basis of Tris-HCl buffers, EDTA and guanidine salts are added, acting as the protein deformers and nuclease inhibitors, making the virus inactive. But this does not affect the integrity of the viral nucleic acid. This facilitates the extraction of nucleic acid and analysis for the subsequent samples, which is safer during inspection and transit but not suitable for isolated cultivation.

Ordering Information

P/N	DESCRIPTION
HX-K60-1	10ml tube with 3ml VTM-N, one oral swab, one nasal swab, one biohazard specimen bag
HX-K60-2	5ml tube with 2ml VTM-N, one oral swab, one nasal swab, one biohazard specimen bag
HX-K60-3	5ml tube with 1ml VTM-N, one oral swab, one nasal swab, one biohazard specimen bag



VTM Collection and Transport Kit






On the basis of Hanks solution, Bovine Serum Albumin V and Virus-stable ingredients such as HEPES are added, maintaining virus activity over a wide temperature range, which facilitates the extraction of nucleic acid for the subsequent samples and the isolated culture of the virus.

- Flocking Swab: $\Phi 2.5 \times 150\text{mm}$ (Stick), 3cm break point for oral swab and 8cm break point for nasal swab
- Transport tube: $\Phi 16 \times 58$ (5ml), $\Phi 16 \times 97 / \Phi 16 \times 101$ (10ml)
- Transport medium: 1ml/tube, 3ml/tube
- Biohazard specimen bag: 4 "x6"

Ordering Information

P/N	DESCRIPTION
HX-K59-1	10ml tube with 3ml VTM, one oral swab, one nasal swab, one biohazard specimen bag
HX-K59-2	5ml tube with 2ml VTM, one oral swab, one nasal swab, one biohazard specimen bag
HX-K59-3	5ml tube with 1ml VTM, one oral swab, one nasal swab, one biohazard specimen bag



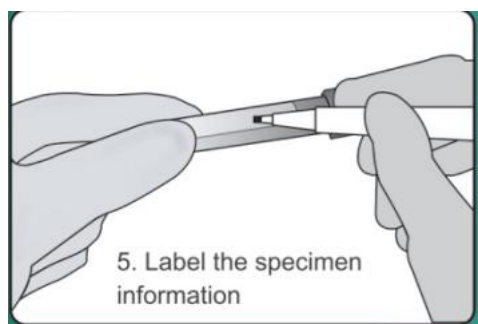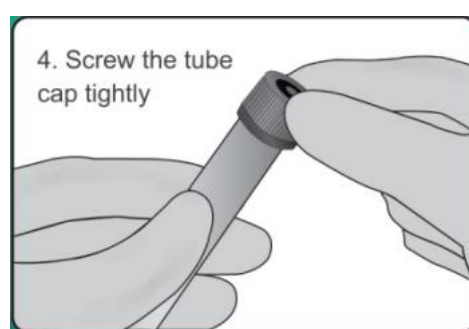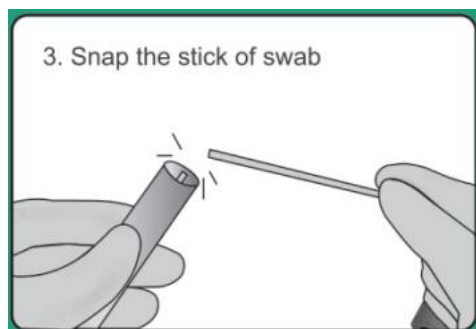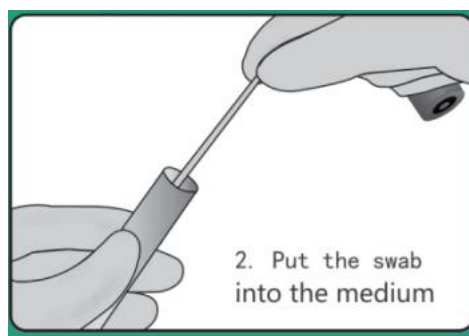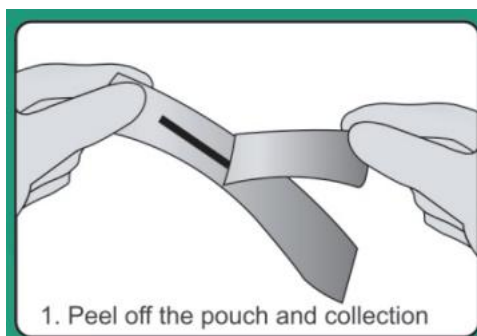

- Disposable flocking swabs, one oral swab, one nasal swab
- VTM and VTM-N Transport media can be chosen as needed.
- Ready to use and easy to tear package, effectively avoid cross contamination.
- Supplied with Biohazard specimen bag, ensure transport safe and reliable.